

Luminex[®]

ARIES[®]

ARIES[®] Systems

ARIES[®] HSV 1&2 Assay (RUO)




The ARIES[®] HSV 1&2 Assay is a rapid method for the direct detection and differentiation of Herpes Simplex Virus 1 and 2 (HSV 1 and HSV 2) DNA in cutaneous or mucocutaneous lesion samples using a sample to answer platform. The ARIES[®] HSV 1&2 RUO Assay offers:

- **Fast Time to Results** — Get answers in less than two hours with minimal hands-on time
- **Flexibility** — Run as many as 12 cassettes per batch
- **Full Integration** — Automate all aspects of testing, from sample extraction through analysis, to enable easy implementation across multiple shifts
- **Error-reducing Safeguards** — Match samples to cassettes with internal barcode scanning, which can help reduce input errors.



ARIES[®] increases research laboratory efficiency, ensures result accuracy, and fits seamlessly into today's lean laboratory.

Workflow

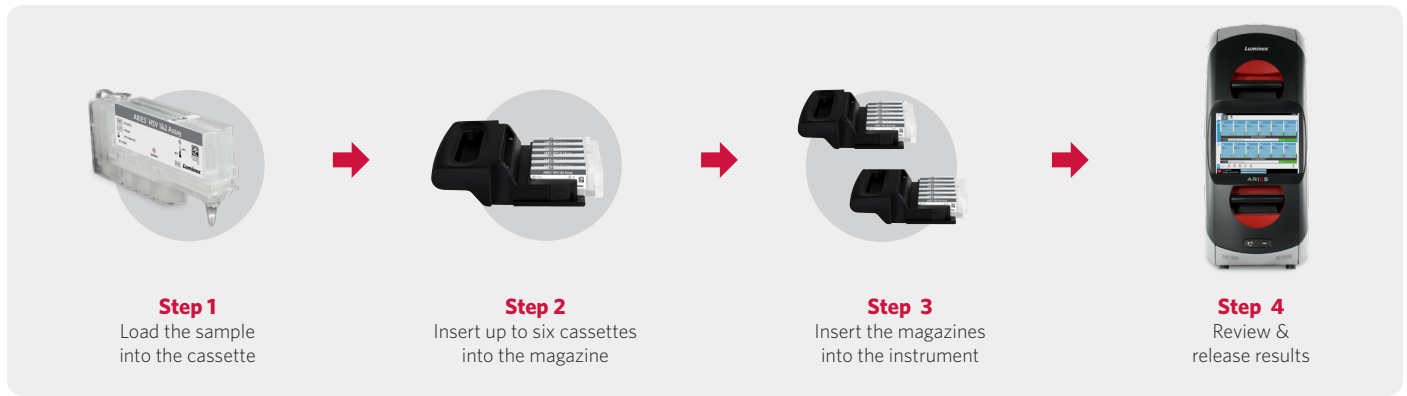

The operator simply adds sample to the sample chamber, puts the cassette in the magazine, loads the magazine into the ARIES® System, and the run will initiate on its own if the auto-run feature has been enabled.

Ordering Information

Description	Part Number
ARIES® HSV 1&2 Assay (RUO) Includes 24 ARIES HSV 1&2 Assay Cassettes	50-10011 (24 tests)
ARIES® HSV 1&2 (RUO) Assay Protocol File Kit	CN-0329-01 (one time order only)
ARIES® Two Module System (RUO) Includes: Instrument System Operation Manual Two Magazines Quick Guide (RUO)	ARIES-M12V1-RUO Two Sample Prep Trays Power Cord Handheld Barcode Scanner
ARIES® M1 System (RUO) Includes: Instrument System Operation Manual One Magazine	ARIES-M6V1-RUO Quick Guide Two Sample Prep Trays Power Cord Handheld Barcode Scanner and Stand

Luminex
complexity simplified.

support@luminexcorp.com or orders@luminexcorp.com

This product is for Research Use Only. Not for use in diagnostic procedures. Please contact Luminex at support@luminexcorp.com to obtain the appropriate product information for your country of residence. The ARIES® System is a class 1(I) laser product. Validation of LIS compatibility must be performed by the end user.

©2019 Luminex Corporation. All rights reserved. Luminex and ARIES are trademarks of Luminex Corporation, registered in the U.S. and other countries.

luminexcorp.com/hsv-ruo/

HEADQUARTERS

UNITED STATES

+1.512.219.8020

info@luminexcorp.com

EUROPE

+31.73.800.1900

europe@luminexcorp.com

CANADA

+1.416.593.4323

info@luminexcorp.com

CHINA

+86.21.8036.9888

info@luminexcorp.com

JAPAN

+81.3.5545.7440

info@luminexcorp.com